

**ORDERING INFORMATION:**

2085 Mobile RPW unit

**Standard Accessories:**

- 28699 - D-ACTOR® Hand piece applicator
- 17638 - R15 15mm ESWT transmitter
- 28724 - D20-S D-ACTOR® 20mm transmitter
- 4248 - Conductor™ transmission gel 250 ml (8.5 oz) bottle
- 28701 - Revision kit that includes the following:
  - 28737 Sealing set
  - 28738 Guide tube (x2)
  - 28736 Projectile (x2)

**Optional Accessories:**

- 28745 - Mobile RPW carrying bag
- 28725 - D35 D-ACTOR® 35 mm transmitter
- 28726 - F15 Focus-Lens 15 mm transmitter
- 28728 - DI15 Deep Impact® 15 mm transmitter
- 28729 - C15 CERAmA-x® ESWT 15 mm transmitter

DJO International  
1a Guildford Business Park  
Guildford  
Surrey  
GU2 8XG  
ENGLAND  
Tel: +44 (0)1483 459 659  
Fax: +44 (0)1483 459 470  
[www.DJOglobal.eu](http://www.DJOglobal.eu)



*Together in Motion™*

## Mobile RPW Shockwave



## MOBILE RPW

Chattanooga's Radial Pressure Wave technology is now mobile. The Mobile RPW is compact and light enough to be easily carried while also delivering a host of features and benefits.

## USER PROFILE

This Radial Pressure Wave clinical unit is ideally suited to the following applications:

- Orthopaedics
- Rheumatology
- Physical Medicine
- Physiotherapy
- Sport Clubs
- Acupuncture
- Aesthetic
- Veterinary
- Podiatrists
- Chiropractors

## INDICATIONS

- Myofascial Trigger Points (MTrP)
  - Localising and deactivating focal points of tenderness
- Activation of muscle and connective tissue
  - Increasing circulation
  - Increasing histological changes
- Insertional tendinopathies (sub chronic - chronic)
  - Plantar Fasciitis, heel spurs
  - Achillodynia
  - Infra-patella tendonosis
  - Chronic Hamstring tendonosis
  - Chronic Hip Trochanteric tendonosis
  - Iliotibial band syndrome
  - Tendinosis calcarea and supraspinatus tendon
  - Medial and Lateral humeral epicondylitis




The Mobile RPW D-Actor Applicator allows quick and efficient uses of applicator transmitters for desired therapeutic effects during treatment.


### D-ACTOR® Applicator


#### Applications:

- Myofascial trigger points
- Other disorders of tendon insertions
- Activation of muscle and connective tissue
- Acupuncture shock wave therapy
- Cellulite treatment (stages I – III)
- Muscle and connective tissue tightening
- Elimination of lymphatic congestion
- Body shaping
- Regulation of muscle tone
- Medical massage

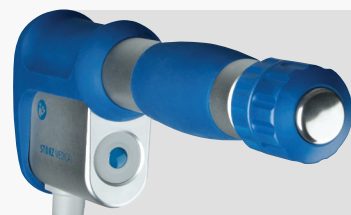

### D-ACTOR® Transmitters

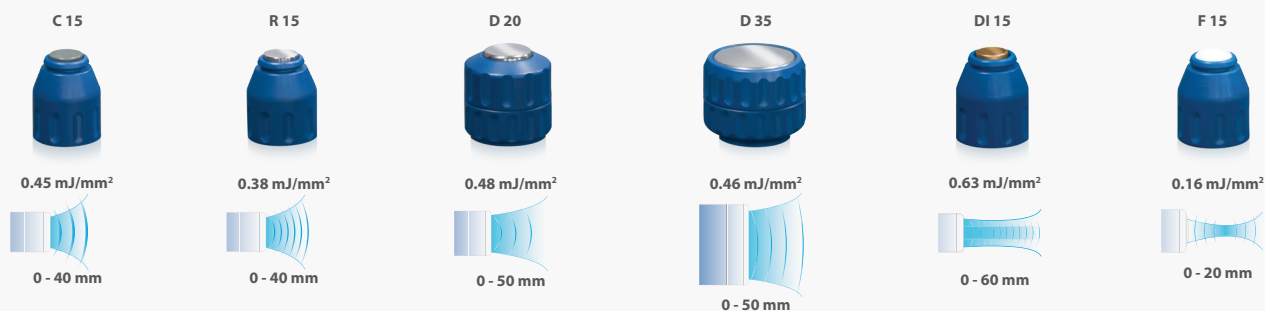

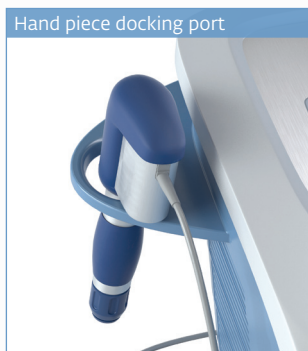

When the need for a powerful yet portable Radial Pressure Wave device is called upon in the clinic or the athletic sports field, the Chattanooga Mobile RPW is the ideal choice. Its lightweight and compact design, together with its quiet, built-in compressor makes it perfect for portable use. User friendliness is at the heart of the Mobile RPW's design. The simple, ergonomically laid out controls enable all of the key parameters to be easily changed, while an illuminated LED display clearly indicates pressure, frequency and number of pulses. Initiation of the pulses is via a trigger mechanism on the hand piece.

#### FACTS AND FIGURES

- Weight of the unit: Approx. 9 kg
- Dimensions (L x W x H): 32 x 32 x 21 cm
- Built in compressor
- Pulse frequency: 1-15 Hz
- Application pressure: 1-4 Bar, 15 MPa
- Oscillating D-Actor



#### USER INTERFACE FEATURES

- |                           |                        |
|---------------------------|------------------------|
| 1 Pressure screen 1-4 Bar | 4 Intensity adjustment |
| 2 Pulse counter screen    | 5 Pulse counter reset  |
| 3 Frequency screen        | 6 Frequency adjustment |